



PROFESSIONAL FITTING AND INFORMATION GUIDE

NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear

CAUTION: US FEDERAL LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A LICENSED PRACTITIONER

TABLE OF CONTENTS

INTRODUCTION	3
DESCRIPTION	3
LENS PARAMETERS	3
ENHANCED SPHERE AND ASPHERE	3
MULTIFOCAL	3
ENHANCED MULTIFOCAL	3
TRANSMITTANCE CURVES	4
ACTIONS	4
INDICATIONS	5
SPHERE AND ASPHERE	5
MULTIFOCAL	5
CONTRAINDICATIONS, WARNINGS, PRECAUTIONS AND ADVERSE REACTIONS	5
GENERAL FITTING GUIDELINES	5
SELECTION OF PATIENTS	5
FITTING GUIDELINES (SPHERE AND ASPHERE)	6
FITTING PROCEDURE	6
CLINICAL ASSESSMENT	7
FITTING GUIDELINES MULTIFOCAL	7
FITTING PROCEDURE	7
INITIAL LENS SELECTION	8
CLINICAL ASSESSMENT	9
TRIAL LENS FITTING	9
ADAPTATION	10
OTHER SUGGESTIONS	10
EVALUATION OF THE LENSES	11
OVER-REFRACTION AND LENS CHANGES	11
IMPROVING VISION	12

INSTRUCTIONS FOR THE MULTIFOCAL WEARER	12
MONOVISION FITTING GUIDELINES	13
SELECTION OF PATIENTS	13
EYE SELECTION	13
SPECIAL FITTING CONSIDERATIONS	14
NEAR ADD DETERMINATION	14
TRIAL LENS FITTING.....	14
ADAPTATION	14
OTHER SUGGESTIONS	15
INSTRUCTIONS FOR THE MONOVISION WEARER.....	15
RECOMMENDED WEARING SCHEDULE.....	16
HANDLING OF LENSES.....	16
LENS PLACEMENT.....	17
LENS REMOVAL	17
PATIENT LENS CARE DIRECTIONS	17
LUBRICATING/REWETTING LENSES ON-EYE	17
CARE FOR A NON-MOVING LENS	17
REPORTING OF ADVERSE REACTIONS.....	18
HOW SUPPLIED.....	18
CAUTION	19
IMPORTANT	19
DESCRIPTION.....	19
TRANSMITTANCE CURVE	19
ACTIONS	20
INDICATIONS	20
CONTRAINDICATIONS (REASONS NOT TO USE)	20
WARNINGS	20
PRECAUTIONS	21
ADVERSE REACTIONS (PROBLEMS AND WHAT TO DO).....	22
FITTING	23
WEARING SCHEDULE.....	23
REPORTING OF ADVERSE REACTION	23
LENS CARE DIRECTIONS	23
HOW SUPPLIED.....	23

INTRODUCTION

NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear are made from etafilcon A with a water content of 58% by weight.

For a complete listing of lens parameters, please refer to LENS PARAMETERS section below.

This Fitting and Information Guide has been developed to provide the Eye Care Professional with information covering characteristics of the NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear, and to illustrate fitting procedures. Please read carefully and keep this information for future use.

See Package Insert for Actions, Contraindications, Warnings, Precautions, Adverse Reactions and Patient Lens Care Directions.

DESCRIPTION

NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear are made of a hydrophilic copolymer of 2-hydroxyethyl methacrylate (2-HEMA) and Methacrylic Acid (MAA). When hydrated, the lenses consist of 42% (etafilcon A) and 58% water by weight when immersed in normal buffered saline. The lens polymer contains a UV absorbing compound with a blue visibility-handling tint, which is color additive “reactive blue19” per 21 CFR part 73.3121. A benzotriazole UV absorbing monomer is used to block UV radiation. The UV Blocking averages 98% in the UVB range of 280 nm to 315 nm and 84% in the UVA range of 316 nm to 380 nm. The etafilcon A name has been adopted by the United States Adopted Names Council (USAN).

WARNING:

UV absorbing contact lenses aren't substitutes for protective UV absorbing eyewear, for example UV absorbing goggles or sunglasses, because they don't completely cover the eye and surrounding area. Patients should continue to use UV absorbing eyewear as directed.

LENS PARAMETERS

NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear are available in sphere, asphere, and multifocal designs, which cover the cornea and a portion of the adjacent sclera with the following dimensions:

Enhanced Sphere and Asphere

- * Diameter: 12.0mm to 15.0mm
- * Center Thickness: 0.08 @ -3.00 D (varies with power)
- * Base Curve: 7.80mm to 10.00mm
- * Powers: +20.00 D to -20.00 D

Multifocal

- * Diameter: 12.0mm to 15.0mm
- * Center Thickness: 0.12 @ -3.00 D (varies with power)
- * Base Curve: 7.80mm to 10.00mm
- * Powers: +20.00 D to -20.00 D
- * ADD Powers: Single universal add power effective up to +3.00 D ADD requirement.

Enhanced Multifocal

- * Diameter: 12.0mm to 15.0mm
- * Center Thickness: 0.08 @ -3.00 D (varies with power)
- * Base Curve: 7.80mm to 10.00mm

* Powers: +20.00 D to -20.00 D

* ADD Powers: Single universal add power effective up to +3.00 D ADD requirement.

The physical/optical properties of the lens are:

Refractive index: 1.402


Light Transmittance: 95 % minimum.

Water Content: 58%

Oxygen Permeability (Dk): $19.73 \times 10^{-11} \text{ (cm}^2/\text{s) x [ml O}_2/\text{ (ml x mmHg)]}$

TRANSMITTANCE CURVES

Typical Transmittance Profile of -3.00 D NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses with UV blocker versus a human cornea from a 24-year-old person and a human crystalline lens from a 25-year-old person.

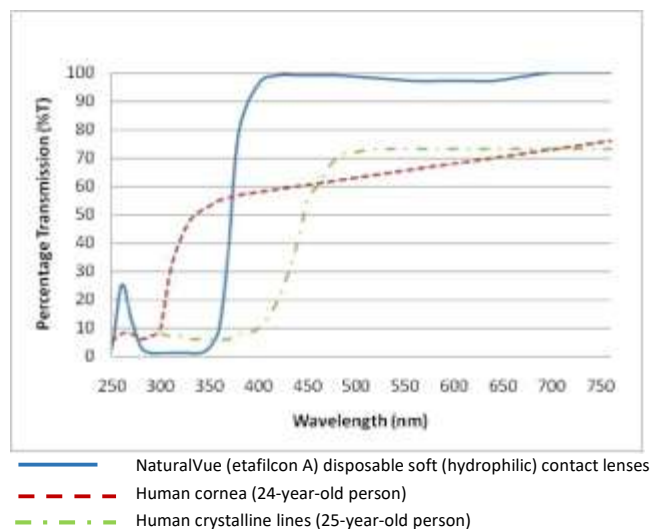

Notes:

1. Lerman, S., Radiant Energy and the Eye, MacMillan, New York, 1980, p. 58, Figure 2-21.

2. Waxler, M., Hitchins, V.M., Optical Radiation and Visual Health, CRC Press, Boca Raton, Florida, 1986, p. 19, Figure 5.

WARNING:

UV absorbing contact lenses aren't substitutes for protective UV absorbing eyewear, for example UV absorbing goggles or sunglasses, because they don't completely cover the eye and surrounding area. You should continue to use UV absorbing eyewear as directed.

ACTIONS

In the hydrated state, the NaturalVue (etafilcon A) Contact Lenses, when placed on the cornea, acts as a refracting medium to focus light rays on the retina.

The visibility tinting of the NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses allows the lens to become visible to the wearer when the lens is not on the eye. The NaturalVue (etafilcon A) Contact Lenses block 84% of UVA radiation and 98% UVB radiation average across the spectrum. (Please refer to accompanying transmittance curve graph)

Note: Long term exposure to UV radiation is a part of the risk factors associated with cataracts. Exposure is according to several factors, for instance environmental conditions (altitude, geography, cloud cover) and personal factors (extent and nature of outdoor activities). UV-absorbing contact lenses help provide protection against harmful UV radiation. However,

clinical studies have not been done to demonstrate that wearing UV-absorbing contact lenses reduces the risk of developing cataracts or other eye disorders. Consult your Eye Care Professional for more information.

INDICATIONS

Sphere and Asphere

NaturalVue (etafilcon A) Enhanced Sphere and Asphere Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear are indicated for daily wear for the correction of ametropia (myopia and hyperopia) in aphakic and/or non-aphakic persons with non-diseased eyes in powers from +20.00 to -20.00 diopters. The powers currently available are +4.00 to -12.00 D. The lenses may be worn by persons who exhibit astigmatism of 2.00 diopters or less that does not interfere with visual acuity.

Multifocal

NaturalVue (etafilcon A) Enhanced Multifocal Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear are indicated for daily wear for the correction of refractive ametropia (myopia and hyperopia), and/or presbyopia in aphakic and/or non-aphakic persons with non-diseased eyes in powers from -20.00 to +20.00 diopters and with non-diseased eyes who may require a reading addition of up to +3.00 D. The powers currently available are +4.00 to -12.25 D. The lenses may be worn by persons who exhibit astigmatism of 2.00 diopters or less that does not interfere with visual acuity.

NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear help protect against transmission of harmful UV radiation to the cornea and into the eye.

The lenses are intended for single-use disposable wear.

CONTRAINDICATIONS, WARNINGS, PRECAUTIONS AND ADVERSE REACTIONS

Refer to Package Insert for CONTRAINDICATIONS, WARNINGS, PRECAUTIONS and ADVERSE REACTIONS.

GENERAL FITTING GUIDELINES

Please see the appropriate sections of this booklet that contain guidelines for Sphere, Asphere, Multifocal, and Monovision fitting techniques.

SELECTION OF PATIENTS

- NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear are indicated for daily wear for vision correction of refractive ametropia (myopia and hyperopia) and presbyopia in aphakic and/or non-aphakic persons with non-diseased eyes.
- The lens is intended for single-use disposable wear.
- Persons who require vision correction and who would not or could not adhere to a recommended wear and care regimen for NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear or are unable to place and remove the lenses should not be provided with them. Failure to follow handling and wearing instructions could lead to serious eye infections that might result in corneal ulcers.
- Patient communication is vital because it relates not only to patient selection, but also to ensuring patient compliance. It is also necessary to discuss the information contained in the Patient Information Booklet with the patient at the time of the initial examination.
- Patients selected to wear NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear should be chosen for their motivation to wear contact lenses, general

health and cooperation. The Eye Care Professional must take care in selecting, examining and instructing contact lens patients. Patient hygiene and willingness to follow practitioner instructions are essential to their success.

- A detailed history is crucial to determining patient needs and expectations. Your patient should be questioned regarding vocation, desired lens-wearing time (full or part-time), and desired lens usage (reading, recreation or hobbies).
- Initial evaluation of the trial lens should be preceded by a complete eye examination, including visual acuity with and without correction at both distance and near, keratometry and slit lamp examination.
- It is normal for the patient to experience mild symptoms such as lens awareness, variable vision, occasional tearing (watery eyes) and slight redness during the adaptation period. Although the adaptation period varies for each individual, generally within one week these symptoms will disappear. If these symptoms persist, the patient should be instructed to contact their Eye Care Professional.

FITTING GUIDELINES (Sphere and Asphere)

FITTING PROCEDURE

- Perform a preliminary evaluation to determine distance refraction as well as to rule out contraindications to contact lens wear as described in the Package Insert.
- Lens power is determined from the patient's spherical equivalent prescription corrected to the corneal plane (i.e., corrected for the vertex distance used for the refraction)
- Make the initial base curve selection if more than one is available. In clinical tests, the NaturalVue (etafilcon A) SPHERE and ASPHERE Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear in the 8.5mm/14.2mm parameters have performed successfully on eyes with a range of keratometry readings. However, corneal curvature measurements should be performed to establish the patient's baseline ocular status. Trial lenses should be placed on each of the patient's eyes and evaluated after the patient has adjusted to the lenses.
- Place the lens on the eye. Allow the lens to remain on the eye long enough to achieve a state of equilibrium. Small variations in the tonicity, pH or the lens solutions and individual tear composition may cause slight changes in fitting characteristics.
- If the initial lens selection covers the patient's cornea fully, provides discernible movement (0.10mm to 0.30mm after blink) or is mobile using the push-up test, is comfortable for the patient and provides satisfactory visual performance, it is a well-fitted lens and can be dispensed. (See Criteria for a Well-Fitted Lens under CLINICAL ASSESSMENT).
- Full coverage of the cornea is defined as the lens edge extending beyond the limbal area in all directions. Initial lens evaluation should be done after at least 10 minutes of lens wear to allow the lens to stabilize and any tearing to subside.
- Following a blink, the lens should move vertically on the patient's eye about 0.10mm to 0.30 mm. Using a slit lamp, this movement can be estimated by comparing it with the one-millimeter lens peripheral bevel width.
- The Push-up Test is also a reliable indicator of a well-fit lens. With the patient looking straight ahead, use your index finger on the patient's lower lid to nudge the edge of the lens upward while observing lens movement, then pull the lid back down and observe the return of the lens. A well-fitted lens will move freely upward, stopping shortly after passing the limbus and then return freely to its original position.

- When lenses are dispensed for vision correction, the wearer must be supplied with an appropriate wearing regimen and must fully understand all lens handling and emergency lens care instructions to prevent lens damage as described in the Package Insert.

CLINICAL ASSESSMENT

- Criteria for a Well-Fitted Lens
 - sufficient lens movement to allow tear exchange under the lens during a blink in either primary or upward gaze.
 - good centration in primary gaze and full corneal coverage in all fields of gaze.
 - satisfactory results on the push-up test.
 - good comfort and a stable visual response (with over-refraction if needed).

- Characteristics of a Tight Lens

A tight lens fit should not be dispensed. A tight lens fit would display some or all of the following characteristics:

- insufficient or no lens movement during a blink in either primary or upgaze.
- unsatisfactory results on the push-up test showing a lens that resists movement, remains decentered or returns slowly to its original position.
- good comfort.
- fluctuating vision between blinks.

- Characteristics of a Loose Lens

A loose lens fit should not be dispensed. A loose lens fit would display some or all of the following characteristics:

- Excessive movement in either primary or upgaze.
- poor centration in primary and upgaze.
- unsatisfactory results on the push-up test showing a lens that moves very easily, remains decentered or returns very quickly to its original position or even drops down lower than its original position.
- reduced comfort.
- vision may be blurred after the blink.

FITTING GUIDELINES Multifocal

FITTING PROCEDURE

- NaturalVue (etafilcon A) Multifocal Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear uses a unique extended depth of focus design, which allows a single universal 'ADD' power to provide optimal intermediate and near vision for a large range of reading additions, while maintaining clear distance vision. NaturalVue (etafilcon A) Multifocal Contact Lenses are easy to fit, as the starting lens powers are the same as the patient's distance spherical contact lens powers (i.e. vertex corrected spherical equivalent lens powers).
- Perform a preliminary evaluation to determine distance and near refraction as well as to rule out contraindications to contact lens wear as described in the Package Insert.
- Ideal Candidates:
 - refractive cylinder ≤ 1.00 D.
 - Emphasis on tasks where it is advantageous to have objects simultaneously in focus over a large range of viewing distances.
 - Expectations consistent with actual everyday visual demands.

- Motivated to wear lenses and understands that vision may not always be as sharp as with spectacles (since spectacles would also generally contain a cylindrical correction) for some distances or lighting conditions.
- Unable to adapt to monovision correction.
- Less than Ideal Candidates:
 - Critical or very fine visual demands at both distance and near.
 - Refractive cylinder > 1.00 D (any axis) in one or both eyes or against-the-rule refractive cylinder > 1.00 D in one or both eyes.
 - Monocular distance acuities poorer than 20/20 with spherical equivalent refractive correction.
 - Myopic anisometropia where the refractive error for one of the two eyes is low (< -1.50 D) and has not been habitually corrected.
 - Pupil size larger (> 5 mm) or smaller (< 3 mm) than the norm for the presbyopic population under natural illumination conditions.
 - Abnormal binocular sensory function (e.g., amblyopia or strabismus).
 - Highly satisfied, adapted monovision or highly satisfied, adapted current multifocal contact lens wearers.
 - Any other contraindications to successful contact lens wear such as tear abnormality or lid margin disease.

Initial Lens Selection

- Start with a current refraction performed at the time of the initial lens fitting:
- ***Best Success Tip:** Ensure the refraction is the sphero-cylindrical distance power that provides optimal vision. DO NOT PUSH PLUS.
- ***Best Success Tip:** To double check that the patient is not over-minused, which is somewhat common in presbyopes, isolate the 20/40 line on a Snellen chart. Binocularly add $+0.25$ D over the distance refraction. This should not blur the 20/40 line noticeably. Add another $+0.25$ D (a total of $+0.50$ D over the distance refraction) and the 20/20 row should just begin to be not quite as clear. Add another $+0.25$ D (a total of $+0.75$ D over the distance refraction), and the 20/40 line should be noticeably blurry. Then use the Duochrome (Red-Green), binocularly in a darkened room, and reduce the plus until the green side letters are first seen as clearer. The final refraction should be one “click” (-0.25 D) into the green.
- Determine ocular dominance using your preferred method. A sensory method is beneficial. The dominance determination will help with any small power adjustments that may be needed.
- The initial diagnostic trial lens power is determined from the patient’s “one click into green” spherical equivalent DISTANCE prescription corrected to the corneal plane (i.e., corrected for the vertex distance used for the refraction). Use of the [NaturalVue QuickStart Calculator](#) is helpful in determining the initial diagnostic lens power. It is important to get the distance vision clear with this lens design.
- Make the initial base curve selection if more than one is available. In clinical tests, NaturalVue (etafilcon A) Multifocal Contact Lenses in the 8.3mm/14.5mm parameters have performed successfully on eyes with a range of keratometry readings. However, corneal curvature measurements should be performed to establish the patient’s baseline ocular status. Trial lenses should be placed on each of the patient’s eyes and evaluated after the patient has adjusted to the lenses.
- Place the lens on the eye. Allow the lens to remain on the eye long enough to achieve a state of equilibrium (generally 10 – 15 minutes). Small variations in the tonicity, pH or the

lens solutions and individual tear composition may cause slight changes in fitting characteristics.

- If the initial lens selection covers the patient's cornea fully, provides discernible movement (0.10mm to 0.30mm after blink) or is mobile using the push-up test, is comfortable for the patient and provides satisfactory visual performance, it is a well fitted lens and can be dispensed. (See Criteria for a Well-Fitted Lens under CLINICAL ASSESSMENT).
- Full coverage of the cornea is defined as the lens edge extending beyond the limbal area in all directions. Initial lens evaluation should be done after at least 10 minutes of lens wear to allow the lens to stabilize and any tearing to subside.
- Following a blink, the lens should move vertically on the patient's eye about 0.10mm to 0.30 mm. Using a slit lamp, this movement can be estimated by comparing it with the one-millimeter lens peripheral bevel width.
- The Push-up Test is also a reliable indicator of a well-fit lens. With the patient looking straight ahead, use your index finger on the patient's lower lid to nudge the edge of the lens upward while observing lens movement, then pull the lid back down and observe the return of the lens. A well-fitted lens will move freely upward, stopping shortly after passing the limbus and then return freely to its original position.
- When lenses are dispensed for vision correction, the wearer must be supplied with an appropriate wearing regimen and must fully understand all lens handling and emergency lens care instructions to prevent lens damage as described in the Package Insert.

CLINICAL ASSESSMENT

- For a good prognosis the patient should have adequately corrected distance and near visual acuity in each eye. The amblyopic patient or the patient with significant astigmatism (greater than 1.00 Diopter) in one eye may not be a good candidate for multifocal correction with the NaturalVue (etafilcon A) Multifocal Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear.
- Occupational and environmental visual demands should be considered. If the patient requires critical vision (visual acuity and stereopsis), it should be determined by trial whether this patient can function adequately with multifocal correction. Multifocal contact lens wear may not be optimal for such activities as:
 - visually demanding situations such as operating potentially dangerous machinery or performing other potentially hazardous activities.
 - driving automobiles (for example, driving at night). Patients who cannot pass their state drivers' license requirements with multifocal correction should be advised to not drive with this correction OR may require that additional over-correction be prescribed.

Trial Lens Fitting

A trial fitting is performed in the office to allow the patient to experience multifocal correction. Lenses are fitted according to the directions in the fitting guidelines.

Case history and standard clinical evaluation procedures should be used to determine the prognosis. With trial lenses of the proper power in place observe the reaction to this mode of correction.

Immediately after the correct power lenses are in place, walk across the room and have the patient look at you. Assess the patient's reaction to distance vision under these circumstances. Then have the patient look at familiar near objects such as a watch face or fingernails.

Again, assess the reaction. As the patient continues to look around the room at both near and

distance objects, observe the reactions. Only after these vision tasks are completed should the patient be asked to read print. Evaluate the patient's reaction to large print (e.g. typewritten copy) at first and then graduate to newsprint and finally smaller type sizes.

After the patient's performance under the above conditions are completed, tests of visual acuity and reading ability under conditions of moderately dim illumination should be attempted. An initial unfavorable response in the office, while indicative of a guarded prognosis, should not immediately rule out a more extensive trial under the usual conditions in which a patient functions.

Adaptation

Visually demanding situations should be avoided during the initial wearing period. A patient may at first experience some mild blurred vision, dizziness, headaches, and a feeling of slight imbalance. You should explain the adaptational symptoms to the patient. These symptoms may last for a brief minute or for several weeks. The longer these symptoms persist, the poorer the prognosis for successful adaptation.

To help in the adaptation process the patient can be advised to first use the lenses in a comfortable familiar environment such as in the home.

Some patients feel that automobile driving performance may not be optimal during the adaptation process. This is particularly true when driving at night. Before driving a motor vehicle, it may be recommended that the patient be a passenger first to make sure that their vision is satisfactory for operating an automobile. During the first several weeks of wear (when adaptation is occurring), it may be advisable for the patient to only drive during optimal driving conditions. After adaptation and success with these activities, the patient should be able to drive under other conditions with caution.

Other Suggestions

The success of the multifocal correction may be further improved by having your patient follow the suggestions below.

- Having supplemental spectacles to wear over the multifocal contact lenses for specific visual tasks may improve the success of the multifocal correction. This is particularly applicable for those patients who cannot meet state licensing requirements with a multifocal correction.
- Make use of proper illumination when carrying out visual tasks.

Success in fitting multifocal corrections can be improved by the following suggestions.

- Refine the lens powers if there is trouble with adaptation. Accurate lens power is critical for presbyopic patients.
- Emphasize the benefits of the clear near vision in straight ahead and upward gaze with multifocal correction.

- **Criteria for a Well-Fitted Lens**

- sufficient lens movement to allow tear exchange under the lens during a blink in either primary or upward gaze.
- good centration in primary gaze and full corneal coverage in all fields of gaze.
- satisfactory results on the push-up test.
- good comfort and a stable visual response (with over-refraction if needed).

- **Characteristics of a Tight Lens**

A tight lens fit should not be dispensed. A tight lens fit would display some or all of the following characteristics:

- insufficient or no lens movement during a blink in either primary or upgaze.

- unsatisfactory results on the push-up test showing a lens that resists movement, remains decentered or returns slowly to its original position.
 - good comfort.
 - fluctuating vision between blinks.
- Characteristics of a Loose Lens
- A loose lens fit should not be dispensed. A loose lens fit would display some or all of the following characteristics:
- Excessive movement in either primary or upgaze.
 - poor centration in primary and upgaze.
 - unsatisfactory results on the push-up test showing a lens that moves very easily, remains decentered or returns very quickly to its original position or even drops down lower than its original position.
 - reduced comfort.
 - vision may be blurred after the blink.

Evaluation of the Lenses

- At the end of the settling period, ask the patient to assess vision at distance, intermediate and near under binocular conditions.
- ***Best Success Tip:** Use real world assessments (e.g., looking out a window, checking cell phone, reading a magazine, reading a watch, etc.).
- Record monocular and binocular distance and near visual acuities.
- Assess lens fit and movement.
- ***Best Success Tip:** The lens is designed to provide a well-centered, minimal movement fit. The lens performs best with < 0.30 mm of decentration and lens movement less than 0.30 mm of either sag or lag. A gentle nudge of the lens with the push-up test should show the lens to move upon push-up even if it is not moving with a blink.

Over-Refraction and Lens Changes

- With this lens design over refraction may not always show a response. It is better to make lens power changes based on patient assessments of vision than on over refraction.
- ***Best Success Tip:** The most important power selection needed to get a successful fit with this design, is the power to obtain clear distance vision. If the distance vision is not clear, add -0.25 D in the dominant eye to clear up distance vision, and another -0.25 D in the dominant eye if needed. If they are still having distance vision complaints, then add -0.25 D in the non-dominant eye. Once distance vision is clear, the intermediate and near vision will improve as well.
- If a lens change is required based on patient assessment that distance vision is not clear, change the lens or lenses and repeat the vision assessment. If vision is acceptable for both distance and near, check lens fit, and proceed with dispensing.
- Vision generally improves over the first 1-2 days with the NaturalVue (etafilcon A) Multifocal lens. A follow-up visit at 2-5 days is recommended to check that the patient is able to perform most daily visual tasks.

Improving Vision

It is important to remember two things with this multifocal contact lens design:

1. The initial diagnostic trial lens power is determined from the patient's "one click into green" spherical equivalent DISTANCE prescription corrected to the corneal plane (i.e., corrected for the vertex distance used for the refraction).
2. The most important power to get a successful fit with this lens design is the power needed to obtain clear distance vision. It is important to get the distance vision clear with this lens design before fine-tuning the intermediate or near vision.
 - Distance Vision and Distance Symptom Complaints: Add -0.25 D, starting with the DOMINANT eye, add another -0.25 D to the DOMINANT eye if needed, and then add -0.25 D to the other eye if needed to further improve distance vision.
 - Near Vision and Near Symptom Complaints: Add +0.25 D, starting with the NON-DOMINANT eye, then another +0.25 D to the NON-DOMINANT eye, and then add +0.25 D to the other eye if needed to further improve near vision.
 - Best Success Tip: Whether one or both lenses are changed, always determine which option provides the best balance between distance and near needs for each individual patient.

The decision to fit a patient with a multifocal correction is most appropriately left to the Eye Care Professional in conjunction with the patient after carefully considering the patient's needs.

All patients should be supplied with a copy of the **Patient Information Booklet**.

INSTRUCTIONS FOR THE MULTIFOCAL WEARER

- You should be aware that as with any type of lens correction, there are advantages and compromises to multifocal contact lens correction. The benefit of clear near vision in straight ahead and upward gaze that is available with multifocal contact lenses may be accompanied by a vision compromise that may reduce your visual acuity and depth perception for distance and near tasks. Some patients have experienced difficulty adapting to it.

Symptoms, such as mild blurred vision, dizziness, headaches and a feeling of slight imbalance, may last for a brief minute or for several weeks as adaptation takes place. The longer these symptoms persist, the poorer your prognosis for successful adaptation. You should avoid visually demanding situations during the initial adaptation period. It is recommended that you first wear these contact lenses in familiar situations, which are not visually demanding. For example, it might be better to be a passenger rather than a driver of an automobile during the first few days of lens wear. It is recommended that you only drive with multifocal correction if you pass your state drivers' license requirements with multifocal correction.

- Some multifocal patients will never be fully comfortable functioning under low levels of illumination, such as driving at night. If this happens, you may want to discuss with your Eye Care Professional having additional contact lenses prescribed so that both eyes are corrected for distance when sharp distance binocular vision is required.
- Some multifocal patients require supplemental spectacles to wear over the multifocal correction to provide the clearest vision for critical tasks. You should discuss this with your Eye Care Professional.

- It is important that you follow your Eye Care Professional's suggestions for adaptation to multifocal contact lenses. You should discuss any concerns that you may have during and after the adaptation period.
- The decision to be fit with multifocal correction is most appropriately left to the Eye Care Professional in conjunction with you, after carefully considering and discussing your needs.

MONOVISION FITTING GUIDELINES

SELECTION OF PATIENTS

- For a good prognosis the patient should have adequately corrected distance and near visual acuity in each eye. The amblyopic patient or the patient with significant astigmatism (greater than 1.00 Diopter) in one eye may not be a good candidate for monovision with the NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear.
- Occupational and environmental visual demands should be considered. If the patient requires critical vision (visual acuity and stereopsis), it should be determined by trial whether this patient can function adequately with monovision. Monovision contact lens wear may not be optimal for such activities as:
 - visually demanding situations such as operating potentially dangerous machinery or performing other potentially hazardous activities.
 - driving automobiles (for example, driving at night). Patients who cannot pass their state drivers' license requirements with monovision correction should be advised to not drive with this correction OR may require that additional over-correction be prescribed.

Eye Selection

Generally, the non-dominant eye is corrected for near vision. The following test for eye dominance can be used.

- Ocular Preference Determination Methods

Method 1: Determine which eye is the "sighting dominant eye." Have the patient point to an object at the far end of the room. Cover one eye. If the patient is still pointing directly at the object, the eye being used is the dominant (sighting) eye.

Method 2: Determine which eye will accept the added power for near with the least reduction in distance vision. Place a trial spectacle near add lens in front of one eye and then the other while the distance refractive error correction is in place for both eyes.

Determine whether the patient functions best with the near add lens over the right or left eye.

- Refractive Error Method

For anisometropic corrections, it is generally best to fit the more hyperopic (less myopic) eye for distance and the more myopic (less hyperopic) eye for near.

- Visual Demands Method

Consider the patient's occupation during the eye selection process to determine the critical vision requirements. If a patient's gaze for near tasks is usually in one direction correct the eye on that side for near.

Example: A secretary who places copy to the left side of the desk will usually function best with the near lens on the left eye.

Special Fitting Considerations

Unilateral Lens Correction

There are circumstances where only one contact lens is required. As an example, an emmetropic patient would only require a near lens while a bilateral myope may require only a distance lens.

Examples:

- A presbyopic emmetropic patient who requires a +1.75 diopter add would have a +1.75 lens on the near eye and the other eye left without a lens.
- A presbyopic patient requiring a +1.50 diopter add who is -2.50 diopters myopic in the right eye and -1.50 diopters myopic in the left eye may have the right eye corrected for distance and the left uncorrected for near.

Near Add Determination

Always prescribe the lens power for the near eye that provides optimal near acuity at the midpoint of the patient's habitual reading distance. However, when more than one power provides optimal reading performance, prescribe the least plus (most minus) of the powers.

Trial Lens Fitting

A trial fitting is performed in the office to allow the patient to experience monovision correction. Lenses are fitted according to the directions in the general fitting guidelines.

Case history and standard clinical evaluation procedure should be used to determine the prognosis. Determine which eye is to be corrected for distance and which eye is to be corrected for near. With trial lenses of the proper power in place observe the reaction to this mode of correction.

Immediately after the correct power lenses are in place, walk across the room and have the patient look at you. Assess the patient's reaction to distance vision under these circumstances. Then have the patient look at familiar near objects such as a watch face or fingernails.

Again, assess the reaction. As the patient continues to look around the room at both near and distance objects, observe the reactions. Only after these vision tasks are completed should the patient be asked to read print. Evaluate the patient's reaction to large print (e.g. typewritten copy) at first and then graduate to newsprint and finally smaller type sizes.

After the patient's performance under the above conditions are completed, tests of visual acuity and reading ability under conditions of moderately dim illumination should be attempted. An initial unfavorable response in the office, while indicative of a guarded prognosis, should not immediately rule out a more extensive trial under the usual conditions in which a patient functions.

Adaptation

Visually demanding situations should be avoided during the initial wearing period. A patient may at first experience some mild blurred vision, dizziness, headaches, and a feeling of slight imbalance. You should explain the adaptational symptoms to the patient. These symptoms may last for a brief minute or for several weeks. The longer these symptoms persist, the poorer the prognosis for successful adaptation.

To help in the adaptation process the patient can be advised to first use the lenses in a comfortable familiar environment such as in the home.

Some patients feel that automobile driving performance may not be optimal during the adaptation process. This is particularly true when driving at night. Before driving a motor vehicle, it may be recommended that the patient be a passenger first to make sure that their

vision is satisfactory for operating an automobile. During the first several weeks of wear (when adaptation is occurring), it may be advisable for the patient to only drive during optimal driving conditions. After adaptation and success with these activities, the patient should be able to drive under other conditions with caution.

Other Suggestions

The success of the monovision technique may be further improved by having your patient follow the suggestions below.

- Having a third contact lens (distance power) to use when critical distance viewing is needed.
- Having a third contact lens (near power) to use when critical near viewing is needed.
- Having supplemental spectacles to wear over the monovision contact lenses for specific visual tasks may improve the success of monovision correction. This is particularly applicable for those patients who cannot meet state licensing requirements with a monovision correction.
- Make use of proper illumination when carrying out visual tasks.

Success in fitting monovision can be improved by the following suggestions.

- Reverse the distance and near eyes if a patient is having trouble adapting.
- Refine the lens powers if there is trouble with adaptation. Accurate lens power is critical for presbyopic patients.
- Emphasize the benefits of the clear near vision in straight ahead and upward gaze with monovision.

The decision to fit a patient with a monovision correction is most appropriately left to the Eye Care Professional in conjunction with the patient after carefully considering the patient's needs.

All patients should be supplied with a copy of the **Patient Information Booklet**.

INSTRUCTIONS FOR THE MONOVISION WEARER

- You should be aware that as with any type of lens correction, there are advantages and compromises to monovision contact lens correction. The benefit of clear near vision in straight ahead and upward gaze that is available with monovision may be accompanied by a vision compromise that may reduce your visual acuity and depth perception for distance and near tasks. Some patients have experienced difficulty adapting to it.

Symptoms, such as mild blurred vision, dizziness, headaches and a feeling of slight imbalance, may last for a brief minute or for several weeks as adaptation takes place. The longer these symptoms persist, the poorer your prognosis for successful adaptation. You should avoid visually demanding situations during the initial adaptation period. It is recommended that you first wear these contact lenses in familiar situations, which are not visually demanding. For example, it might be better to be a passenger rather than a driver of an automobile during the first few days of lens wear. It is recommended that you only drive with monovision correction if you pass your state drivers' license requirements with monovision correction.

- Some monovision patients will never be fully comfortable functioning under low levels of illumination, such as driving at night. If this happens, you may want to discuss with your Eye Care Professional having additional contact lenses prescribed so that both eyes are corrected for distance when sharp distance binocular vision is required.

- Some monovision patients require supplemental spectacles to wear over the monovision correction to provide the clearest vision for critical tasks. You should discuss this with your Eye Care Professional.
- It is important that you follow your Eye Care Professional's suggestions for adaptation to monovision contact lens wear. You should discuss any concerns that you may have during and after the adaptation period.
- The decision to be fit with monovision correction is most appropriately left to the Eye Care Professional in conjunction with you, after carefully considering and discussing your needs.

RECOMMENDED WEARING SCHEDULE

The prescribing Eye Care Professional for each individual patient, based upon a full examination and patient history as well as the practitioner's experience and professional judgment should determine the wearing schedule. Patients should be given a wearing schedule and carefully instructed on the handling and care of their lenses as discussed in the Package Insert. Also be sure to complete the personal wearing/replacement schedule record in the patient information booklet. Eye Care Professionals may prescribe the lens for single-use daily disposable wear. The lens is intended for single-use, daily disposable wear. (See the factors discussed in the WARNINGS section.)

Follow up examinations are necessary to ensure continued successful contact lens wear and to ascertain the effects of the lenses on the eyes. The following schedule is a suggested guideline for daily wear contact lenses:

- 24 hours post-dispensing
- 7 days
- 1 month
- 3 months
- every 6 months thereafter

HANDLING OF LENSES

- Always wash, rinse and dry hands before handling contact lenses.
- All traces of soap, perfumes, hair sprays, creams and lotions should be removed from your hands and around the eyes.
- Keep your nails trimmed and clean – long fingernails can tear or split contact lenses.
- Avoid picking up lenses with your fingernails – they can only be held safely between the fingertips or with soft plastic tipped tweezers specially designed for contact lenses.
- Do not touch your lenses with sharp, pointed objects (fingernails, pens, etc.).
- Make sure lenses do not get caught on the edge of blister pack.
- Never allow lenses to dry out and never try to insert them when they are in a dehydrated state.
- Do not use the lens if the pack is open or damaged.
- Never use expired lenses or solutions.
- **Verify that the lens is right side out.** The lens should assume a natural, curved, bowl-like shape. If the lens edges tend to point outward, the lens is inside out. Another method is to gently squeeze the lens between the thumb and forefinger. The edges should turn inward. If the lens is inside out, the edges will turn slightly outward.

Water Activity

Instruction for Use:

- Do not expose your contact lenses to water while you are wearing them.

WARNING:

Water can harbor microorganisms that can lead to severe infection, vision loss or blindness. If your lenses have been submersed in water when swimming in pools, lakes or oceans, you should discard them and replace them with a new pair. Ask your eyecare practitioner (professional) for recommendations about wearing your lenses during any activity involving water.

LENS PLACEMENT

To avoid confusion, make it a habit of inserting the same lens first. Check the lens is clean and moist – if not, rinse with saline or multipurpose solution. Check the lens is not inside out. And then follow these steps to apply the lens to your eye:

- Balance the lens on the tip of your forefinger.
- Pull the lower lid down with the middle finger of the same hand.
- Hold the upper lid firmly from above with the middle finger of the other hand. Make sure the finger is placed just where the eyelid meets the eyelashes.
- The eye is now wide open and the cornea exposed.
- Place the lens directly on the eye. Move your eyes around to help center the lens.
- Release the lower and then the upper lid. Look down and gently close the eyes. Rub the top lid lightly to rub out any trapped air bubbles.
- If the lens is not properly centered on the cornea, gently manipulate by using fingertips and eyelids.

LENS REMOVAL

- Rewet the contact lenses with a few drops of solution to help the lens slide off the eye, especially if they are a little dry.
- Place index finger on the lens.
- Look up.
- Slide lens to white part of the eye using index finger.
- Gently squeeze the lens between thumb and forefinger to remove the lens.

PATIENT LENS CARE DIRECTIONS

Refer to Package Insert for Lens Care Directions.

LUBRICATING/REWETTING LENSES ON-EYE

The Eye Care Professional may recommend a lubricating/rewetting solution that can be used to wet (lubricate) lenses while they are being worn to make them more comfortable.

CARE FOR A NON-MOVING LENS

Refer to Package Insert for Care for A Sticking (Nonmoving) Lens.

REPORTING OF ADVERSE REACTIONS

If a patient experiences any serious adverse events associated with the wear of NaturalVue (etafilcon A) Daily Soft (Hydrophilic) Contact Lenses, please notify:

Visioneering Technologies, Inc.

30 Mansell Court, Suite 215

Alpharetta, GA 30076

1-844-VTI-LENS

(1-844-884-5367)

HOW SUPPLIED

Each lens is supplied sterile in a blister pack containing buffered saline solution with Tween 80, Sodium Hyaluronate, and Polyethylene Glycol. The blister is labeled with the base curve, diameter, diopter power, manufacturing lot number, and expiration date of the lens. (ADD, cylinder and axis will be included as appropriate.)



NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses for Daily Wear

CAUTION

US Federal law restricts this device to sale by or on the order of a licensed practitioner.

IMPORTANT

This package Insert is intended for the Eye Care Professional but should be made available to patients upon request. The Eye Care Professional should provide the patient with appropriate instructions that pertain to the patient's prescribed lenses and recommended wearing schedule.

DESCRIPTION

The NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses are available in a spheric, aspheric, and multifocal designs.

NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses are made of a hydrophilic copolymer of 2-hydroxyethyl methacrylate (2-HEMA) and Methacrylic Acid (MAA). When hydrated, the lenses consist of 42% (etafilcon A) and 58% water by weight when immersed in normal buffered saline. The lens polymer contains a UV absorbing compound with a blue visibility-handling tint, color additive "reactive Blue19," 21 CFR part 73.3121. A benzotriazole UV absorbing monomer is used to block UV radiation. The UV Blocking averages 98% in the UVB range of 280 nm to 315 nm and 84% in the UVA range of 316 nm to 380 nm. The etafilcon A name has been adopted by the United States Adopted Names Council (USAN).

LENS PROPERTIES

Refractive index: 1.402

Light Transmittance: 95% minimum

Water Content: 58%

Oxygen Permeability (Dk):

$$19.73 \times 10^{-11} \text{ (cm}^2/\text{s) x [ml O}_2\text{/ (ml x mmHg)]}$$

LENS PARAMETERS

NaturalVue Enhanced Sphere and Asphere

* Diameter: 12.0mm to 15.0mm

* Center Thickness: 0.08 @ -3.00 D (varies with power)

* Base Curve: 7.80mm to 10.00mm

* Powers: +20.00 D to -20.00 D

NaturalVue Multifocal

* Diameter: 12.0mm to 15.0mm

* Center Thickness: 0.120 @ -3.00 D (varies with power)

* Base Curve: 7.80mm to 10.00mm

* Powers: +20.00 D to -20.00 D

* ADD Powers: Single universal add power effective up to +3.00 D

NaturalVue Enhanced Multifocal

* Diameter: 12.0mm to 15.0mm

* Center Thickness: 0.08 @ -3.00 D (varies with power)


* Base Curve: 7.80mm to 10.00mm

* Powers: +20.00 D to -20.00 D

* ADD Powers: Single universal add power effective up to +3.00 D

TRANSMITTANCE CURVE

Typical Transmittance Profile of -3.00 D NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses with UV blocker versus a human cornea from a 24-year-old person and a human crystalline lens from a 25-year-old person.

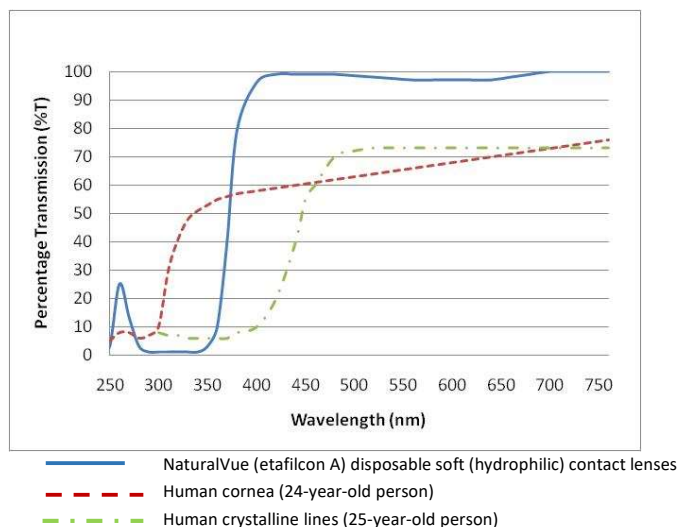

Notes:

1. Lerman, S., Radiant Energy and the Eye, MacMillan, New York, 1980, p. 58, Figure 2-21.
2. Waxler, M., Hitchins, V.M., Optical Radiation and Visual Health, CRC Press, Boca Raton, Florida, 1986, p. 19, Figure 5.

WARNING:

UV absorbing contact lenses aren't substitutes for protective UV absorbing eyewear, for example UV absorbing goggles or sunglasses, because they don't completely cover the eye and surrounding area. You should continue to use UV absorbing eyewear as directed.

ACTIONS

In the hydrated state, the NaturalVue (etafilcon A) Contact Lenses, when placed on the cornea, act as a refracting medium to focus light rays on the retina. The visibility tinting of the NaturalVue (etafilcon A) Contact Lenses allow the lens to become visible to the wearer when the lens is not on the eye. The NaturalVue (etafilcon A) Contact Lenses block 84% of UVA radiation and 98% UVB radiation average across the spectrum. (Please refer to accompanying transmittance curve graph)

Note: Long term exposure to UV radiation is a part of the risk factors associated with cataracts. Exposure is according to a number of factors, for instance environmental conditions (altitude, geography, cloud cover) and personal factors (extent and nature of outdoor activities). UV-absorbing contact lenses help provide protection against harmful UV radiation. However, clinical studies have not been done to demonstrate that wearing UV-absorbing contact lenses reduces the risk of developing cataracts or other eye disorders. Consult your Eye Care Professional for more information.

INDICATIONS

Sphere and Asphere

NaturalVue (etafilcon A) Sphere and Asphere Daily Disposable Soft (Hydrophilic) Contact Lenses are indicated for daily wear for the correction of ametropia (myopia and hyperopia) in aphakic and/or non-aphakic persons with non-diseased eyes in powers from +20.00 to -20.00 diopters. The powers currently available are +4.00 D to -12.00 D. The lenses may be worn by persons who exhibit astigmatism of 2.00 diopters or less that does not interfere with visual acuity.

Multifocal

NaturalVue (etafilcon A) Multifocal Daily Disposable Soft (Hydrophilic) Contact Lenses are indicated for daily wear for the correction of refractive ametropia (myopia and hyperopia), and/or presbyopia in aphakic and/or non-aphakic persons with non-diseased eyes in powers from -20.00 to +20.00 diopters and with non-diseased eyes who may require a reading addition of up to +3.00 D. The powers currently available are +4.00 D to -12.25 D. The lenses may be worn by persons who exhibit astigmatism of 2.00 diopters or less that does not interfere with visual acuity.

NaturalVue (etafilcon A) Daily Disposable Soft (Hydrophilic) Contact Lenses help protect against transmission of harmful UV radiation to the cornea and into the eye.

The lenses are intended for single-use disposable wear.

CONTRAINDICATIONS (REASONS NOT TO USE)

DO NOT USE the NaturalVue Lens when any of the following conditions exist:

- Acute and subacute inflammation or infection of the anterior chamber of the eye.
- Any eye disease, injury, or abnormality that affects the cornea, conjunctiva, or eyelids.
- Severe insufficiency of lacrimal secretion (dry eyes).
- Corneal hypoesthesia (reduced corneal sensitivity).
- Any systemic disease that may affect the eye or be exaggerated by wearing contact lenses.
- Allergic reactions of ocular surfaces or adnexa (surrounding tissue) that may be induced or exaggerated by wearing contact lenses or use of contact lens solutions.
- Any active corneal infection (bacterial, fungal, or viral).
- If eyes become red or irritated.
- Use of any medication that is contraindicated or interferes with contact lens wear, including ocular medications.
- Allergy to any ingredient in the contact lens material.
- History of recurrent eye or eyelid infections, adverse effects associated with contact lens wear, intolerance or an unusual response to contact lens wear.
- The patient is unable or unwilling to follow the Eye Care Professionals' directions for removal and disposal of the lenses or unable to obtain assistance to do so.

WARNINGS

Advise patient of the following warnings pertaining to contact lens wear:

- Problems with contact lenses could result in serious injury to the eye. It is essential that patients follow their Eye Care Professional's directions and all labeling instructions for proper use of lenses.
- Eye problems, including a sore or lesion on the cornea (corneal ulcers), can develop rapidly and lead to loss of vision.
- Studies have shown that contact lens wearers who are smokers have a higher risk incidence of adverse reactions than nonsmokers, especially when lenses are worn overnight or while sleeping.
- If a patient experiences eye discomfort, such as foreign body sensation, excessive tearing, vision changes, or redness of the eye or other problems, the patient should immediately remove lenses and promptly contact his or her Eye Care Professional.
- Daily wear lenses are not indicated for overnight wear, and patients should be instructed not to wear lenses while sleeping. Clinical study results have shown that the risk of serious adverse reactions is increased when lenses are worn overnight.
- Non-compliance with the manufacture's labeled lens care instruction may put the patient at significant risk of developing a serious eye infection.

- Tap water, distilled water, homemade saline solutions or saliva should NOT be used at any time with contact lenses. The use of tap and distilled water has been associated with Acanthamoeba keratitis, a corneal infection that is resistant to treatment and cure.

Water Activity

Instruction for Use:

- Do not expose your contact lenses to water while you are wearing them.

WARNING:

Water can harbor microorganisms that can lead to severe infection, vision loss or blindness. If your lenses have been submersed in water when swimming in pools, lakes or oceans, you should discard them and replace them with a new pair. Ask your eyecare practitioner (professional) for recommendations about wearing your lenses during any activity involving water.

PRECAUTIONS

Special Precautions for Eye Care Professional:

- Due to the small number of patients enrolled in clinical investigation of lenses, all refractive powers, design configurations, or lens parameters available in the lens material are not evaluated in significant numbers. Consequently, when selecting an appropriate lens design and parameters, the Eye Care Professional should consider all characteristics of the lens that can affect lens performance and ocular health, including oxygen permeability, wettability, central, and peripheral thickness, and optic zone diameter.
- The potential impact of these factors on the patient's ocular health should be carefully weighed against the patient's need for refractive correction; therefore, the continuing ocular health of the patient and lens performance on eye should be carefully monitored by the prescribing Eye Care Professional.
- Fluorescein, a yellow dye, should not be used while the patient is wearing the lenses, because the lenses will absorb this dye and become discolored. Whenever fluorescein is used in eye, flush the eyes with sterile saline solution. Wait at least 10 minutes before reinserting the lenses. If it is not possible to flush the eyes, wait at least 1 hour before wearing the lenses. If inserted too soon, the lenses may absorb remaining fluorescein.
- Before leaving the Eye Care Professional's office, the patient should be able to promptly remove lenses or

should have somebody else available who can remove the lenses for him or her.

- Eye Care Professionals should instruct the patient to remove the lenses immediately if the eyes become red or irritated.

Eye Care Professionals should carefully instruct patients about the following safety precautions, including the need for routine eye examinations being necessary to help assure the continued health of the patient's eyes:

- If the lens sticks (stops moving) on the eye, follow the recommended directions in Care for a Sticking Lens. The lens should move freely on the eye for the continued health of the eye. If non-movement of the lens continues, you should immediately consult your Eye Care Professional.
- Always wash and rinse hands before handling lenses. Do not get cosmetics, lotions, soaps, creams, deodorants, aerosol products or hair sprays in the eyes or on the lenses. It is best to put on lenses before putting on make-up. Water-base cosmetics are less likely to damage lenses than oil-base products.
- Do not touch contact lenses with the fingers or hands if the hands are not free of foreign materials, as microscopic scratches of the lenses may occur, causing distorted vision or injury to the eye.
- Exposure to water while wearing contact lenses in activities such as swimming, water skiing, and hot tubs may increase the risk of ocular infection, including but not limited to Acanthamoeba keratitis.
- Ask the Eye Care Professional about wearing contact lenses during sporting activities.
- Never wear lenses beyond the period recommended by the Eye Care Professional.
- If aerosol products such as hair spray are used while wearing lenses, exercise caution and keep eyes closed until the spray has settled.
- Always handle lenses gently and avoid dropping them.
- Avoid all harmful or irritating vapors and fumes while wearing lenses.
- Aphakic patients should not be fitted with lenses until the determination is made that the eye has healed completely from surgery.
- Inform the Eye Care Professional about being a contact lens wearer.
- Never use tweezers or other tools to remove lenses from the lens blister pack unless specifically indicated for that use. Pour the lens into the hand.
- Do not touch the lens with fingernails.
- Always discard lenses after the recommended wearing schedule prescribed by the Eye Care Professional.
- Always contact the Eye Care Professional before using any medicine in the eyes.
- Always inform the employer of being a contact lens wearer. Some jobs may require use of eye protection equipment or may require that the patient not wear contact lenses.

- As with any contact lens, follow-up visits are necessary to ensure the continuing health of the patient's eyes. The patient should be instructed as to a recommended follow-up schedule.
- Patients should never exceed the prescribed wearing schedule regardless of how comfortable the lenses feel. Doing so increases the risk of adverse effects.
- Do not use lenses past the expiration date.
- Certain medications may cause dryness of the eye, increased lens awareness, lens intolerance, blurred vision or visual changes. These include, but are not limited to, antihistamines, decongestants, diuretics, muscle relaxants, tranquilizers, and those for motion sickness. Caution patients using such medications accordingly and prescribe proper remedial measures.
- Patients who are pregnant or oral contraceptive users could develop visual change or change in lens tolerance when using lenses.
- Do not use if the sterile blister package is opened or damaged.
- Diabetics may have reduced corneal sensitivity and may be more prone to corneal injury and do not heal as quickly or completely as non-diabetics.
- Patients who wear contact lenses to correct presbyopia may not achieve the best correct visual acuity for either far or near vision. Vision requirements vary with the individual and should be considered when selecting the most appropriate type of lens for each patient.
- Patients should be instructed to never allow anyone to wear their lenses. They have been prescribed to fit their eyes and to correct their vision to the degree necessary. Sharing lenses greatly increase the chance of eye infections.

ADVERSE REACTIONS (Problems and What To Do)

The patient should be informed that the following problems may occur when wearing contact lenses:

- Your eye stinging, burning, itching (irritation), or other eye pain.
- Comfort is less compared to when lens was first placed on eye.
- There may be a feeling of something in the eye (foreign body, scratched area).
- Excessive watering (tearing) of the eyes.
- Unusual eye secretions.
- Redness of the eyes.
- Reduced sharpness of vision (poor visual acuity).
- Blurred vision, rainbows, or halos around objects.
- Sensitivity to light (photophobia).
- Feeling of dryness.
- Foreign body sensation.

If the patient notices any of the above, he or she should be instructed to:

- Immediately remove lenses.
- If the discomfort or problem stops, then look closely at the lens. If the lens is in any way damaged, do not

put the lens back on the eye. If the lens has dirt, an eyelash, or other foreign body on it, or the problem stops and the lens appears undamaged; the patient should thoroughly clean and rinse the lens, then reinsert it. After reinsertion, if the problem continues, the patient should immediately remove the lens and consult the Eye Care Professional.

- If the above symptoms continue after removal of the lens, or upon insertion of a new lens, the patient should immediately remove the lens and contact his or her Eye Care Professional, who must determine the need for examination, treatment or referral without delay. (See Important Treatment Information for Adverse Reactions.) A serious condition such as infection, corneal ulcer, corneal vascularization, or iritis may be present, and may progress rapidly. Less serious reactions such as abrasions, epithelial staining or bacterial conjunctivitis must be managed and treated carefully to avoid more serious eye damage. Additionally, contact lens wear may be associated with ocular changes which require consideration of discontinuation or restriction of wear. These include but are not limited to local or generalized corneal edema, epithelial microcysts, epithelial staining, infiltrates, neovascularization, endothelial polymegathism, tarsal papillary changes, conjunctival injection or iritis.

- Patients should be reminded to keep a spare pair of lenses with them at all times.

Important Treatment Information for Adverse Reactions

Sight-threatening ocular complications associated with contact lens wear can develop rapidly, and therefore early recognition and treatment of problems are critical. Infectious corneal ulceration is one of the most serious potential complications, and may be ambiguous in its early stage. Signs and symptoms of infectious corneal ulceration include discomfort, pain, inflammation, purulent discharge, sensitivity to light, cells and flare and corneal infiltrates.

Initial symptoms of a minor abrasion and an early infected ulcer are sometimes similar. Accordingly, such epithelial defect, if not treated properly, may develop into an infected ulcer. In order to prevent serious progression of these conditions, a patient presenting symptoms of abrasions or early ulcers should be evaluated as a potential medical emergency, treated accordingly, and be referred to a corneal specialist when appropriate. Standard therapy for corneal abrasions such as eye patching or the use of steroids or steroid/antibiotic combinations may exacerbate the condition. If the patient is wearing a contact lens on the affected eye when examined, the lens should be removed immediately.

FITTING

- The lens must move adequately on the eye for a suitable fit. If the fit is judged to be too tight, the patient must be refit into a lens which provides the criteria of a well-fitted lens.
- Fitting techniques for other contact lenses may not be applicable to the fitting of these lenses.
- Copies of Professional Fitting and Information Guide for NaturalVue (etafilcon A) Contact Lenses are available without charge at www.vtivision.com or from:
Visioneering Technologies, Inc.
30 Mansell Court, Suite 215
Alpharetta, GA 30076
Tel: 1-844-VTI-LENS x102
(1-844-884-5367) x102
www.vtivision.com

WEARING SCHEDULE

It is recommended that contact lens wearers see their Eye Care Professional twice each year or if directed, more frequently.

Daily Wear (Less than 24 hours, while awake):

There may be a tendency for the **NEW** daily wear patient to over-wear the lenses initially. Initial daily wearing schedule should be stressed to these patients.

The wearing schedule should be determined by the Eye Care Professional. The maximum suggested wearing time is:

DAY	1	2	3	4	5	6
HOURS	6	8	10	12	14	All waking hours

NaturalVue Lenses are only indicated for daily wear (less than 24 hours, while awake). NaturalVue Lenses should NOT be worn for 24 hours or greater, including while asleep, as studies have not been completed to show that the lenses are safe to wear during sleep.

REPORTING OF ADVERSE REACTION

If any serious adverse event or patient experience associated with the wear of NaturalVue (etafilcon A) Daily Soft (Hydrophilic) Contact Lenses occurs, please notify:

Visioneering Technologies, Inc.

30 Mansell Court, Suite 215

Alpharetta, GA 30076

Tel: 1-844-VTI-LENS x102

(1-844-884-5367 x102)

www.vtivision.com

LENS CARE DIRECTIONS

- **Daily disposable** lenses are discarded upon removal from the eye **each day**.
- Patients should be instructed to have a spare pair of lenses with them at all times.

Care for a Sticking (Nonmoving) Lens:

If the lens sticks (cannot be moved), you should use a lubricating or rewetting solution in your eye. You should apply 3 to 4 drops of the recommended lubricating or rewetting solution directly to the eye and wait until the lens begins to move freely on the eye before removing it. If non-movement of the lens continues after 15 minutes, you should IMMEDIATELY consult your Eye Care Professional.

Storage:

- Unopened lenses are sterile and need not be disinfected until the container is broken or the seal has been damaged.

Emergencies:

The patient should be informed that if chemicals of any kind (household products, gardening solutions, laboratory chemicals, etc.) are splashed into the eyes, the patient should:

FLUSH EYES IMMEDIATELY WITH TAP WATER AND THEN REMOVE LENSES PROMPTLY. CONTACT THE EYE CARE PROFESSIONAL OR VISIT A HOSPITAL EMERGENCY ROOM WITHOUT DELAY.

HOW SUPPLIED

Each lens is supplied sterile in a blister pack containing buffered saline solution. Enhanced Sphere and Enhanced Multifocal additionally contain Tween 80, Sodium Hyaluronate, and Polyethylene Glycol in the blister pack buffered saline solution. The blister is labeled with the base curve, diameter, diopter power, manufacturing lot number, and expiration date of the lens.

	Manufactured by: Visioneering Technologies, Inc. 30 Mansell Court, Suite 215 Alpharetta, GA 30076 Tel: 1-844-VTI-LENS x102 1-844-884-5367 x102 www.vtivision.com
---	---

SYMBOL	DESCRIPTION
	Consult instructions for use
	Batch code
	Sterilized Using Steam or Dry Heat
	Use by Date Expressed as: CCYY-MM or CCYY-MM-DD CCYY-MM
	Product to be sold by or on the order of a licensed practitioner.
	UV-Blocking (Indicates the lens is UV blocking)
B.C.	Base Curve (product property)
D	Diopter (Lens Power, product property)
DIA.	Diameter (product property)
	Do not re-use

SYMBOL	DESCRIPTION
	Manufacturer
	Do not use if package is damaged
ADD	ADD (product property)

Revision date January 2024